


Enema bottle

Enema bottle for rectal products includes a polyethylene (PE) return flow regulation valve inserted in the cannula base. This ensures the integrity of the liquid treatment by not allowing it to flow back into the bottle during dispensing.



OVERVIEW

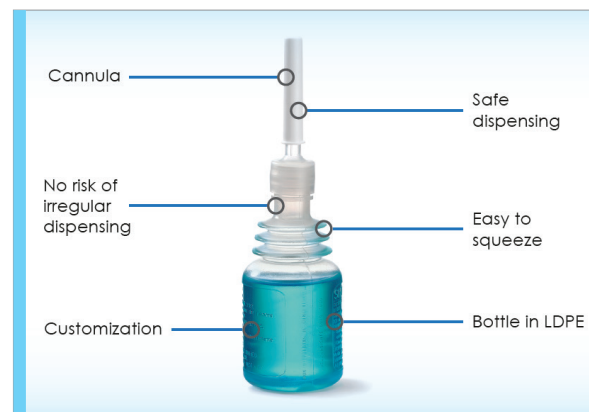

- **Bottle in low density polyethylene (LDPE)** available in 35, 60, 75 and 120ml.



- **Safe dispensing:** The cannula with membrane provides tamper-resistance for the container and contents.



- **Cannula can be pre-assembled with the closure cap.**



- **Silicone treatment of the cannula** upon request.



- **Return flow regulation valve polyethylene (PE)** inserted in the base of the cannula.



- **Customization with embossing and colored caps available.**



- **Translucent bottles allow the fluid level to be visible.**



- **Easy to squeeze:** Soft bottle with thin sidewalls for easy squeezability and total product evacuation.

MARKETS



Pharmaceutical



Rectal



Medical Device



For additional information, please visit:
tekni-plex.com/healthcare